

# BladderScan® BVI 6400

## Mobile Bladder Volume Instrument



### BladderScan®

#### The Standard of Care for Bladder Volume Measurement

The BladderScan® family of products consists of portable ultrasound instruments that measure bladder volume. They provide quick, accurate and reliable results.

Because they are noninvasive and easy to use, BladderScan® instruments are an excellent alternative to diagnostic urinary catheterisation and a perfect tool to help healthcare providers noninvasively diagnose, manage and treat urinary outflow dysfunction.

#### BladderScan® Benefits

- Prevents unnecessary catheterisation
- Helps prevent undiagnosed urinary retention
- Reduces rates of urinary tract infection
- Helps caregivers manage and treat incontinence
- Eliminates unnecessary trauma to patients, preserves their dignity and helps improve patient care
- Cost effective - saves staff time and improves job satisfaction

#### Mobile BladderScan® BVI 6400

The handheld, portable BladderScan® BVI 6400:

- Measures bladder volume noninvasively
- Improves efficiency of health care professionals who examine multiple patients on their rounds
- Performs and stores up to ten exams
- Provides for voice annotation on each exam (10 seconds to record patient ID and relevant exam information) ensuring that valuable patient and exam data is retained
- Allows for exam results and images to be downloaded, viewed and printed using the Verathon™ ScanPoint® imaging service
- Provides fast, accurate and reliable results
- Is easy to use
- Is battery operated and can be calibrated online via ScanPoint®



UROLOGY	ACUTE CARE
<p>For the urologist, BladderScan® applications assist in:</p> <ul style="list-style-type: none"> <li>• Diagnosis of severe bladder outlet obstruction (BOO), lower urinary tract symptoms (LUTS), overactive bladder (OAB), and prostatitis</li> <li>• Monitoring post-operative recovery</li> <li>• Measurement of post-void residual bladder volume (PVR)</li> <li>• Evaluation of urinary retention</li> </ul>	<p>BladderScan® minimises catheterisation and dramatically lowers the incidence of nosocomial urinary tract infection (NUTI) at your facility. Department of Health statistics show that NUTI:</p> <ul style="list-style-type: none"> <li>• Increases the cost of care for hospitals and rehabilitation centres</li> <li>• Lengthens hospital stays by 1- 3 days</li> <li>• Can be life-threatening for patients</li> </ul>
PRIMARY CARE	EXTENDED CARE
<p>BladderScan® aids in the diagnosis and treatment of:</p> <ul style="list-style-type: none"> <li>• Benign prostatic hyperplasia (BPH)</li> <li>• Bladder outlet obstruction (BOO)</li> <li>• Prostatitis</li> <li>• Overactive bladder (OAB)</li> </ul> <p>The BladderScan® BVI 6400 will provide you with the information you need to more easily evaluate urological conditions, so you can prescribe appropriate treatment sooner.</p>	<p>When utilised as part of a comprehensive incontinence management programme, the BladderScan®:</p> <ul style="list-style-type: none"> <li>• Minimizes incontinent episodes</li> <li>• Reduces infection rates (UTI)</li> <li>• Improves patients' quality of life</li> <li>• Saves staff time and increases job satisfaction</li> <li>• Reduces the cost of care</li> </ul>


Noninvasive, Accurate, Reliable, & Easy to Use



# BladderScan® BVI 6400

## Mobile Bladder Volume Instrument

### How BladderScan® Works

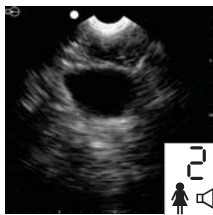

Sample Display - BVI 6400

The BladderScan® calculates bladder volume using patented V<sub>MODE</sub>® technology. V<sub>MODE</sub>® ultrasound is easy to use and comfortable for the patient. When you press the scan button, within seconds, the V<sub>MODE</sub>® technology measures ultrasonic reflections on multiple planes inside the body and produces a 3-dimensional image. Based on this image, the BladderScan® calculates and displays the bladder volume, and displays multiple B-mode images.

Volume measurements made with V<sub>MODE</sub>® ultrasound are more accurate than those from conventional 2-dimensional ultrasound, as they are based on a more complex, 3-D image of the bladder.



By implementing the BladderScan® BVI 6400, you can improve the quality of care for patients with urinary problems.



#### DIAGNOSE

- Measure post-void residual (PVR)
- Differentiate urological problems more efficiently
- Diagnose urinary retention
- Identify blocked Foley catheter
- Verify empty bladder
- Determine hydration status
- Differentiate between types of incontinence to determine appropriate care

#### MANAGE AND TREAT

- Establish volume (vs. time) based intermittent catheterisation protocols
- Discontinue Foley catheter use
- Evaluate need to catheterise after Foley catheter removal
- Establish voiding trials
- Evaluate need to catheterise during intermittent catheterisation
- Establish voiding schedules
- Assist in bladder retraining (biofeedback)

#### PREVENT

- Eliminate unnecessary catheterisation
- Reduce rates of nosocomial UTI
- Minimise antibiotic use associated with UTI
- Minimise catheter-related damage to upper urinary tract
- Prevent undiagnosed (chronic) urinary retention
- Prevent post-operative (acute) urinary retention
- Prevent bladder overdistension
- Reduce incontinent episodes

### Specifications - BVI 6400

Range:	Bladder volume range 0 to 999 ml	Ultrasound Output Parameters:	Maximum SPTA* Intensity: 1.04mW/cm <sup>2</sup> Maximum SPPA* Intensity: 65.0mW/cm <sup>2</sup> Mechanical Index (MI): 0.925 maximum Ultrasound Frequency: 3.7 MHz Scan angle: 120 degrees Mode: V <sub>MODE</sub> ® (multiple, aligned B-mode images)
Accuracy:	The following accuracy specification assumes usage per instructions, scanning a Verathon™ Tissue Equivalent Phantom: 0 to 999 ml ± 15%, ± 15 ml	Operating Conditions:	Temperature: +10° C to +40° C Humidity: 30% to 75%, non-condensing
Scan Time:	Less than 5 seconds	* SPTA = Spatial temporal average; SPPA = Spatial peak pulse average	
Weight:	Less than 11 oz. (309 grams)		
Power:	3.8v Li Ion rechargeable battery		
Display:	Liquid crystal		

The BladderScan® BVI 6400 is CE Marked in accordance with the Medical Device Directive and the Verathon™ quality system is Quality System Certified to ISO 13485:2003 standards. US Patent No. 4926871, 5235985, 6569097, 6884217. INTL Patents Pending. BladderScan®, ScanPoint®, V<sub>MODE</sub>®, Verathon™ and Verathon Medical™ are either registered trademarks or trademarks of Verathon Inc. in the USA and/or other countries. Copyright ©2006. All rights reserved.



**Verathon Medical (United Kingdom) Ltd.**

The Granary Manor Farm Courtyard / Aston Sandford, Aylesbury  
Buckinghamshire, HP17 8JB / United Kingdom  
+44.1844.299.207 / Fax: +44.1844.299.218

[www.verathon.com](http://www.verathon.com)